



Distribuito in ITALIA da
Li StarFish S.r.l.
Via Cavour, 35
20063 Cernusco S/N (MI)
telefono 02-92150794
fax 02-92157285
info@listarfish.it
www.listarfish.it

en

Read entire protocol before use.

1,25(OH)₂ Vitamin D ELISA

I. INTENDED USE

Immunoenzymetric assay for the *in vitro* quantitative measurement of 1,25(OH)₂ Vitamin D in serum.

II. GENERAL INFORMATION

Catalog number : VIT-KAP1921 : 96 tests

III. CLINICAL BACKGROUND

A. Biological activities

Vitamin D is mainly synthesized in the skin from 7-dehydrocholesterol and is partially from dietary and supplementation origin. In the liver, Vitamin D is hydroxylated on carbon 25 to produce the intermediate 25OH Vitamin D. 25OH Vitamin D is further metabolized before it can carry out the functions of Vitamin D on intestine, kidneys, bone and other organs and tissues. This subsequent reaction takes place in the kidneys and in other tissues. Thus 25OH Vitamin D is further hydroxylated in the 1 α -position to produce 1 α ,25-dihydroxyvitamin D (1,25(OH)₂ Vitamin D). In addition to the above-mentioned tissues, placenta of pregnant women and macrophage cells in case of sarcoidosis can also produce some amount of 1,25(OH)₂ Vitamin D.

1,25(OH)₂ Vitamin D is the active form of Vitamin D with regard to the known functions whereas 25OH Vitamin D and Vitamin D itself can be excluded as being physiologically functional. 1,25(OH)₂ Vitamin D stimulates the intestinal absorption of both calcium and phosphorus. It also stimulates bone resorption and mineralization thereby preventing the development of rickets and osteomalacia.

1,25(OH)₂ Vitamin D is also active in other tissues responsible for Calcium transport (placenta, kidney, mammary gland,...) and endocrine glands such as parathyroid glands. 1,25(OH)₂ Vitamin D is rapidly metabolized and its half-life is approximately 12h in plasma. Its main metabolite is calcitroic acid, a C-23 carboxylic derivative, essentially without any biological activity. In addition to this pathway, 1,25(OH)₂ Vitamin D undergoes 24-hydroxylation to produce 1,24,25-trihydroxyvitamin D. This compound has less biological activity than its parent and this metabolic route is considered as a minor pathway.

The levels of 1,25(OH)₂ Vitamin D in plasma or serum is 100 to 1000 less than that of 25OH Vitamin D. Due to its low concentrations and the presence of many similar metabolites, the measurement of 1,25(OH)₂ Vitamin D requires extraction and separation by chromatography.

B. Clinical application

The measurement of circulating 1,25(OH)₂ Vitamin D is indicated in several disorders affecting calcium metabolism such as : phosphate diabetes, sarcoidosis, renal failure, hyper and hypo-parathyroidism, rickets, tumor-associated hypercalcemia, hypercalciuria, Vitamin-resistant dysfunction and treatment with anti-convulsive medication.

IV. PRINCIPLES OF THE METHOD

Only samples and controls, not the calibrators, are extracted with a mixture of solvents and applied on cartridges to separate 1,25(OH)₂ Vitamin D from the other Vitamin D metabolites. After elution of the 1,25(OH)₂ Vitamin D from the samples and controls cartridges, the calibrators, eluted samples and eluted controls are incubated directly in microtiterplate coated with anti-1,25(OH)₂ Vitamin D antibodies.

After an overnight incubation at 4°C, the microtiter plate is washed and the working conjugate solution is then added and incubated for 1 hour at 4°C.

The microtiterplate is then washed to stop the competition reaction. The Chromogenic solution (TMB) is added and incubated for 15 minutes at room temperature. The reaction is stopped with the addition of Stop Solution and the microtiterplate is read at the appropriate wavelength.

The amount of 1,25(OH)₂ Vitamin D is determined colourimetrically by measuring the absorbance, which is inversely proportional to the 1,25(OH)₂ Vitamin D concentration.

A calibration curve is plotted and the 1,25(OH)₂ Vitamin D concentrations of the samples are determined by dose interpolation from the calibration curve.

V. REAGENTS PROVIDED


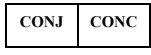
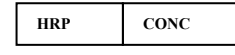
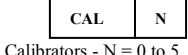
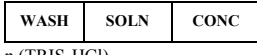
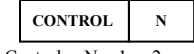
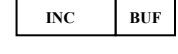
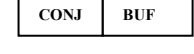

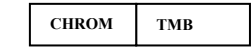
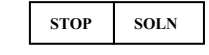
Reagents	96 Test Kit	Colour Code	Reconstitution
 Microtiterplate with 96 breakable wells coated with anti-1,25(OH) ₂ Vitamin D antibodies	96 wells	blue	Ready for use
 1,25(OH) ₂ Vitamin D Concentrated Conjugate	1 vial 1 mL	blue	Dilute 40 x with conjugate buffer
 Concentrated HRP	1 vial 0.2 mL	yellow	Dilute 200 x with conjugate buffer
 Calibrators - N = 0 to 5 (see exact values on vial labels) in phosphate buffer with bovine casein and gentamycin	6 vials lyophilised	yellow	Add 1mL distilled water
 Wash solution (TRIS-HCl)	1 vial 10 mL	brown	Dilute 200 x with distilled water (use a magnetic stirrer)
 Controls - N = 1 or 2 in human plasma with gentamycin	2 vials lyophilised	silver	Add 2mL distilled water
 Incubation Buffer with proclin	1 vial 20 mL	green	Ready for use
 Conjugate Buffer with casein and proclin	1 vial 30 mL	red	Ready for use
 Elution Solution: contains methanol	1 vial 25 mL	white	Ready for use
 Chromogenic solution TMB (Tetramethylbenzidine)	1 vial 25 mL	brown	Ready for use
 Stop solution HCl 1.5 N	1 vial 12 mL		Ready for use

PLATE	COVER		
Adhesive Strips		4	
GEL		42	Store at R.T.
Silica cartridges			

Note : Use Calibrator 0 for dilution of samples with values above the highest calibrator (dilute before extraction step).

VI. SUPPLIES NOT PROVIDED

The following material is required but not provided in the kit:

- 1 Distilled water
- 2 Diisopropylether ("for analysis"; GC purity ≥ 99%)
- 3 Cyclohexane ("for analysis"; GC purity ≥ 99.5 %)
- 4 Ethyl acetate("for analysis"; GC purity ≥ 99.5 %)
- 5 Ethanol absolute ("for analysis"; GC purity ≥ 99.9 %)
- 6 Dichloromethane ("for analysis"; GC purity ≥ 99.8 %)
- 7 Pipettes for delivery of: 50µL, 100µL, 150µL, 200 µL, 1mL and 2 mL (the use of accurate pipettes with disposable plastic tips is recommended)
- 8 Glass tubes (12 x 75 mm) for extraction and for elution (closed with a cap for the extraction step).
- 9 Glass tubes (16 x 100 mm) or (12 x 120 mm), or polypropylene tubes (falcon 2097), for the washing of the cartridges.
- 10 Vortex mixer
- 11 Magnetic stirrer
- 12 Centrifuge operating at 800 g
- 13 Microtiterplate reader capable of reading at 450 nm and 650 (bichromatic reading)

VII. REAGENT PREPARATION

A. **Calibrators:** Reconstitute the calibrators with 1mL distilled water (**just before incubation step**).

B. **Controls:** Reconstitute the controls with 2mL distilled water.

C. **Working HRP conjugate solution:**

! The working HRP conjugate solution is to be prepared absolutely one hour before adding the solution in the plate !

Prepare an adequate volume of working HRP conjugate solution by mixing concentrated conjugate, concentrated HRP and conjugate buffer according to the number of used strips, as indicated in the below table: for example for 6 strips (48 wells): 250µL of concentrated conjugate and 50µL of concentrated HRP to 10mL of conjugate buffer.

Use a vortex to homogenize.

Keep the working HRP conjugate at room temperature and avoid direct sunlight or use a brown glass vial for its preparation.

Nb of strips	Volume of Concentrated Conjugate (µL)	Volume of Concentrated HRP (µL)	Volume of Conjugate Buffer (mL)
1	75	15	3
2	125	25	5
3	150	30	6
4	200	40	8
5	225	45	9
6	250	50	10
7	300	60	12
8	350	70	14
9	400	80	16
10	450	90	18
11	500	100	20
12	550	110	22

D. **Working Wash solution:** Prepare an adequate volume of Working Wash solution by adding 199 volumes of distilled water to 1 volume of Wash Solution (200x). Use a magnetic stirrer to homogenize. Discard unused Working Wash solution at the end of the day.

E. **Extraction solvent:** 2 ml for each control or sample to be tested are needed. **Prepare a fresh solution** of diisopropylether, cyclohexane and ethyl acetate: 50/40/10 volume/volume).

Be careful : exact proportion of volume have to be strictly respected.

- F Washing solvent** : 1 ml for each control or sample to be tested is needed.
Prepare a fresh solution of diisopropylether, cyclohexane, ethyl acetate and ethanol absolute (50/40/10/1 volume/volume).
Be careful : exact proportion of volume have to be strictly respected .

VIII. STORAGE AND EXPIRATION DATING OF REAGENTS

- Before opening or reconstitution, all kits components are stable until the expiry date, indicated on the label, if kept at 2 to 8°C; **except the cartridges which must be stored at room temperature.**
- After their reconstitution, **calibrators** are stable for 4 weeks at 2-8°C. For longer storage periods, aliquots should be made and kept at -20°C for 4 months maximum. Avoid subsequent freeze-thaw cycles.
- After their reconstitution, **controls** are stable for 3 days at 2-8°C. For longer storage periods, aliquots should be made and kept at -20°C for 1 month maximum. Avoid subsequent freeze-thaw cycles.
- Freshly prepared Working Wash solution should be used on the same day.
- After its use, discard working HRP conjugate.
- **Use freshly prepared extraction solvent and washing solvent, do not store them.**
- Alterations in physical appearance of kit reagents may indicate instability or deterioration.

IX. SPECIMEN COLLECTION AND PREPARATION

- The kit is suitable for serum samples.
- Serum samples must be kept at 2-8°C.
- If the test is not run within 24 hrs, storage in aliquots, at -20°C is recommended.
- Avoid subsequent freeze-thaw cycles.
- After thawing, the samples should be vortexed and centrifuged.

X. PROCEDURE

I. Extraction step : ! Only for controls and samples !

1. Label glass tubes (12x75 mm) for extraction: 2 controls and up to 40 samples.
2. Add 0.5mL control or sample in the respective tubes.
3. Dispense 2mL extraction solvent in each tube.
4. Tubes are closed with a cap and placed on a shaker for 1 hour at 1200 rpm.
5. Centrifuge each tube for 5 minutes at room temperature (at 800 g).
6. Supernatants are needed for the next step of separation.

II. Separation step : ! Only for controls and samples !

1. Label glass tubes (16 x 100 mm) or (12 x 120 mm), or polypropylene tubes (falcon 2097), for washing cartridges: 2 controls and up to 40 samples.
2. Put one silica cartridge in each tube.
3. Apply 1.6mL of supernatant (2 x 0.8mL), obtained after extraction step, on cartridge. Let draw by gravity.
4. Wash cartridges with 1mL washing solvent (cf: reagent preparation).
! Be careful: never apply vacuum on cartridges, just let solvent draw by gravity.
5. Add 500µL dichloromethane on each cartridge, let draw by gravity.
6. Add 500µL of distilled water on each cartridge and centrifuge each tube for 5 minutes at room temperature (at 800 g).
7. Label glass tubes (12 x 75 mm) for elution of 1,25(OH)₂ Vitamin D. After centrifugation, transfer cartridges in the corresponding glass tubes.
8. Apply 300µL elution solution on each cartridge to elute 1,25(OH)₂ Vitamin D and centrifuge for 5 minutes at room temperature (at 800 g).
9. **Vortex** the eluted fraction.

Note : After this step, samples must be incubated in coated microtiterplate immediately to avoid degradation.

III. Incubation step :

1. Select the required number of strips for the run. The unused strips should be resealed in the bag with a desiccant and stored at 2-8°C.
2. Secure the strips into the holding frame.
3. Vortex briefly reconstituted calibrators, extracted controls and extracted samples.
4. Pipette 150µL of Incubation Buffer into all wells.
5. Pipette 50µL of each Calibrator (not extracted), eluted controls and eluted samples into the appropriate wells.
6. Incubate for 18±2 hours, at 2-8°C. Cover the plate with a lid or a sealing film.

Prepare the Working HRP conjugate solution 60 min +/- 15 min before washing the wells after the overnight incubation.

7. Aspirate the liquid from each well.
Wash the plate 3 times by:
 - dispensing 0.35mL of Wash Solution into each well and
 - aspirating the content of each well
8. Pipette 200µL of Working HRP conjugate solution into each well.

9. Incubate for 1 hour at 4°C. Cover the plate with a lid or a sealing film.
10. Aspirate the liquid from each well.
Wash the plate 3 times by:
 - dispensing 0.35 ml of Wash Solution into each well and
 - aspirating the content of each well
11. Pipette 200µL of the Chromogenic solution into each well within 15 minutes following the washing step.
12. Incubate the microtiterplate for 15 minutes at room temperature, avoid direct sunlight.
13. Pipette 100µL Stop Solution into each well.
14. Read absorbances at 450 nm (reference filter 630nm or 650 nm) within 1 hour and calculate the results as described in section XI.

XI. CALCULATION OF RESULTS

1. Read the plate at 450 nm against a reference filter set at 650 nm (or 630 nm).
2. Calculate the mean of duplicate determinations.
3. We recommend the use of computer assisted methods to construct the calibration curve. 4-parameter logistic function curve fitting is the preferred method. Reject obvious outliers.
4. By interpolation of the sample OD values, determine the 1,25(OH)₂ Vitamin D concentrations of the samples from the calibration curve.

XII. TYPICAL DATA

The following data are for illustration only and should never be used instead of the real time calibration curve.

1,25(OH) ₂ Vitamin D ELISA	OD units
Calibrator :	
0 pg/mL	2.93
3 pg/mL	2.52
12 pg/mL	1.85
50 pg/mL	1.11
120 pg/mL	0.57
180 pg/mL	0.36

Note : 1 pg/mL = 2.4 pmol/L

XIII. PERFORMANCE AND LIMITATIONS

A. Detection Limit

Twenty zero calibrators were assayed along with a set of other calibrators. The detection limit, defined as the apparent concentration two standard deviations below the average OD at zero binding, was 0.8 pg/mL.

B. Specificity

Cross reactivity of the 1,25(OH)₂ Vitamin D ELISA assay was determined by testing sera with spiked and unspiked cross reactants. The results are summarized in the following table:

Compound and concentration	Cross-Reactivity (%)
1,25(OH) ₂ -Vitamin.D3 at 200 pg/mL	114
1,25(OH) ₂ -Vitamin.D2 at 200 pg/mL	108
25OH-Vitamin-D3 at 1µg /mL	0.004
25OH-Vitamin-D2 at 1µg /mL	0.0003
24,25(OH) ₂ -Vitamin.D3 at 200 ng/mL	0.03
25,26(OH) ₂ -Vitamin.D3 at 400ng/mL	0.02

The effect of potential interfering substances on samples using the DIAsource 1,25(OH)₂ Vitamin D ELISA test was evaluated. Different levels of Hemoglobin, Bilirubin (conjugated and unconjugated), Triglyceride and Vitamin C in serum samples were tested on samples with different 1,25(OH)₂ Vitamin D Concentration. Our acceptance criteria was to have interference of less than 15%. The tested substances did not affect the performance of the DIAsource 1,25(OH)₂ Vitamin D ELISA.

Substance	1,25(OH) ₂ VitaminD (pg/mL)	Concentration of Interferent (mg/dL)	Mean % Variation
Hemoglobin	31.8	250	5.0%
		500	
	186.5	250	
		500	
Bilirubin Conjugated	31.8	50	-12.3%
	186.5	50	
Bilirubin Unconjugated	31.8	50	-0.4%
		100	
	186.5	50	
		100	
Triglyceride	31.8	50	-1.0%
		100	
		250	
	186.5	50	
		100	
		250	
Vitamin C	31.8	100	4.9%
		1000	
	186.5	100	
		1000	

C. Precision

INTRA-ASSAY PRECISION

INTER-ASSAY PRECISION

Sample	N	<X> ± SD (pg/mL)	CV (%)	Sample	N	<X> ± SD (pg/mL)	CV (%)
A	13	18.3 ± 2.6	14.2	A	13	26.7 ± 3.5	13.2
B	13	168.9 ± 8.7	5.1	B	13	83.4 ± 14.6	17.5

SD: Standard Deviation; CV: Coefficient of variation

D. Accuracy

The sample was diluted with Calibrator 0, before extraction step.

DILUTION TEST

Sample Dilution	Theoretical Concentration (pg/mL)	Measured Concentration (pg/mL)	Slope	Y-Intercept	R	Recovery (%)
1/1	118.9	118.9	0.99	1.12	0.99	100
1/2	59.5	60.7				102
1/4	29.7	29.3				99
1/8	14.9	16.6				112

Conversion factor :

From pg/mL to pmol/L : x 2.4

From pmol/L to pg/mL : x 0.42

To the best of our knowledge, no international reference material exists for this parameter.

RECOVERY TEST

Added 1,25(OH) ₂ -Vitamin D (pg/mL)	Recovered 1,25(OH) ₂ -Vitamin D (pg/mL)	Recovered (%)
52.4	54.1	103
104.7	111.1	106
157.1	155.8	99

XIV. INTERNAL QUALITY CONTROL

- If the results obtained for Control 1 and/or Control 2 are not within the range specified on the vial label, the results cannot be used unless a satisfactory explanation for the discrepancy has been given.
- If desirable, each laboratory can make its own pools of control samples, which should be kept frozen in aliquots. Controls which contain azide will interfere with the enzymatic reaction and cannot be used.
- Acceptance criteria for the difference between the duplicate results of the samples should rely on Good Laboratory Practises
- It is recommended that Controls be routinely assayed as unknown samples to measure assay variability. The performance of the assay should be monitored with quality control charts of the controls. It is good practise to check visually the curve fit selected by the computer

XV. REFERENCE INTERVALS

These values are given only for guidance; each laboratory should establish its own normal range of values.

Normal samples tested with DIASource Elisa assay were measured between 19.3 and 53.8 pg/ml.

Patients with renal failure (n = 20) were measured < 6,9 pg/ml.

XVI. PRECAUTIONS AND WARNINGS

Safety

For *in vitro* diagnostic use only.

The human blood components included in this kit have been tested by European approved and/or FDA approved methods and found negative for HBsAg, anti-HCV, anti-HIV-1 and 2. No known method can offer complete assurance that human blood derivatives will not transmit hepatitis, AIDS or other infections. Therefore, handling of reagents, serum or plasma specimens should be in accordance with local safety procedures.

All animal products and derivatives have been collected from healthy animals. Bovine components originate from countries where BSE has not been reported. Nevertheless, components containing animal substances should be treated as potentially infectious.

Avoid any skin contact with all reagents, Stop Solution contains HCl. In case of contact, wash thoroughly with water.

Do not smoke, drink, eat or apply cosmetics in the working area. Do not pipette by mouth. Use protective clothing and disposable gloves.

XVII. BIBLIOGRAPHY

1. Bouillon R.A., Auwerx J.D., Lissens W.D. and Pelemans W.K. (1987). **Vitamin D status in elderly; season substrate deficiency causes 1,25-dihydroxycholecalciferol deficiency.** Am. J. Clin. Nutr., 45:755-763.
2. Iqbal, S.J. (1994). **Vitamin D metabolism and the clinical aspects of measuring metabolites.** Ann. Clin. Biochem., 31:109-124.
3. Mawer E.B. (1980). **Clinical implications of measurements of circulating vitamin D metabolites.** Clinics in Endocr. Metabol., 9:63-79.
4. Jongen M.J.M., Van Ginkel F.C., Vander Vijgh W.J.F., Kuiper S., Netelenbos J.C. and Lips P; (1984). **An international comparison of Vitamin D metabolites measurements.** Clin. Chem., 30:399-403.
5. Deluca H.F. (1979). **The Vitamin D system in the regulation of calcium and phosphorus metabolism.** Nutritional Rev., 37:161-193.
6. Haussler, M.R., McCain, T.A. (1977). **Basic and clinical concepts related to Vitamin D metabolism and action.** N. Engl. J. Med., 297:974-983.
7. Pike J.W., Lee S.M., Meyer M.B. (2014). **Regulation of gene expression by 1,25-dihydroxyvitamin D3 in bone cells: exploiting new approaches and defining new mechanisms.** BoneKEy Reports, 3:482.
8. Hollis B.W. (2010). **Assessment and Interpretation of Circulating 25-Hydroxyvitamin D and 1,25-Dihydroxyvitamin D in the Clinical Environment.** Endocrinol. Metab. Clin. North Am., 39(2): 271–contents.

XVIII.

SUMMARY OF THE PROTOCOL

	CALIBRATORS μL	SAMPLE(S) CONTROLS μL
<u>EXTRACTION</u> Calibrators Samples / Controls Extraction solvent	- - -	- 500 2000
Shaking Centrifugation	1 hour at 1200 rpm 5 minutes at 800 g	
<u>SEPARATION</u> Supernatant from extraction step	-	1600
<u>CARTRIDGE</u> Supernatant Washing Solvent Dichloromethane Distilled water Centrifugation Elution solution Centrifugation	1600μL 1000μL 500μL 500μL 5 minutes at 800 g 300μL 5 minutes at 800 g Vortex	
<u>INCUBATION STEP</u> <u>In microtiterplate</u> Incubation Buffer Calibrators Extracted samples	150μL 50μL -	150μL - 50μL
	Cover the plate with a lid or sealing film Incubate 18 ± 2 h (overnight) at 4°C (2-8°C) Prepare working HRP solution 1 hour before next step Aspirate the contents of each well Wash 3 times with 350μL of Wash Solution and aspirate	
Working HRP Conjugate	200μL	200μL
	Cover the plate with a lid or sealing film and Incubate for 1 hour at 4°C (2-8°C). Aspirate the contents of each well. Wash 3 times with 350 μl of Wash Solution and aspirate.	
TMB	200μL	200μL
	Incubate for 15 min at room temperature (18°C to 25°C).	
Stop Solution	100μL	100μL
	Read on a microtiterplate reader. Record the absorbance of each well at 450 nm (versus 630 or 650 nm).	



Distribuito in ITALIA da
Li StarFish S.r.l.
 Via Cavour, 35
 20063 Cernusco S/N (MI)
 telefono 02-92150794
 fax 02-92157285
 info@listarfish.it
 www.listarfish.it

P.I. Number :
1700551/en

Revision nr :
151020/1

Revision date : 2015-10-20