



Distribuito in ITALIA da
Li StarFish S.r.l.
Via Cavour, 35
20063 Cernusco S/N (MI)
telefono 02-92150794
fax 02-92157285
info@listarfish.it
www.listarfish.it

Bovine IFN- γ ELISA development kit

Product Code: 3119-1A-6

CONTENTS, development kit for 6 plates:

Vial 1 (green top)

Monoclonal antibody MT17.1 (300 μ l)

Concentration: 0.5 mg/ml

Vial 2 (red top)

Biotinylated monoclonal antibody MT307 (50 μ l)

Concentration: 0.5 mg/ml

Vial 3 (white top)

Streptavidin-Alkaline Phosphatase (80 μ l)

Vial 4

Recombinant bovine IFN- γ standard (0.5 μ g)

To ensure total recovery of stated quantity, vials have been overfilled.

STORAGE:

Shipped at ambient temperature. On arrival box 1 should be stored refrigerated at 4-8°C and box 2 should be stored frozen at -20°C.


General

Intended use: For quantitative determination of native and recombinant bovine IFN- γ in solution, e.g. cell culture supernatant and serum/plasma samples. The two mAbs cross-react with native ovine IFN- γ .

Reagents: Antibodies are supplied in sterile-filtered (0.2 μm) PBS with sodium azide (0.02%). Streptavidin-ALP is supplied in 0.1 M Tris buffer with 0.002% Kathon CG.

Standard range: 5-500 pg/ml.

Standard calibration: No international standard exists for calibration.

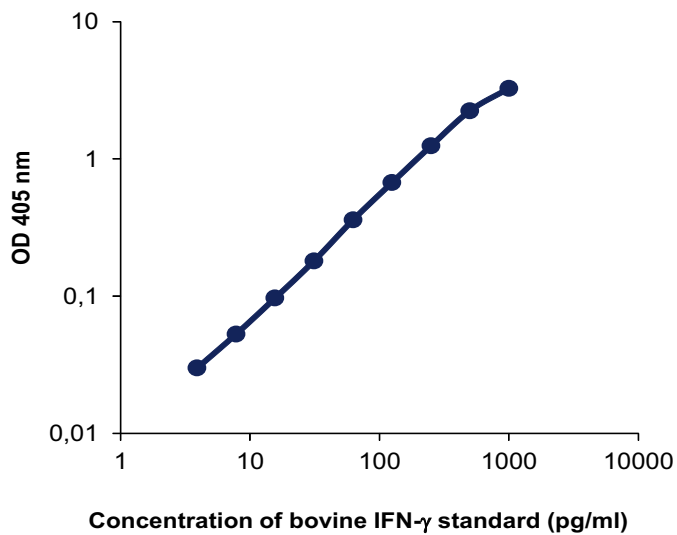

Guidelines for Bovine IFN- γ ELISA

- Day 1**
1. Coat a high protein binding ELISA plate with mAb MT17.1, diluted to 2 $\mu\text{g}/\text{ml}$ in PBS, pH 7.4, by adding 100 $\mu\text{l}/\text{well}$. Incubate overnight at 4-8°C.
- Day 2**
2. Wash twice with PBS (200 $\mu\text{l}/\text{well}$).
 3. Block plate by adding 200 $\mu\text{l}/\text{well}$ of PBS with 0.05% Tween 20 containing 0.1% BSA (incubation buffer). Incubate for 1 hour at room temperature.
 4. Wash five times with PBS containing 0.05% Tween.
 5. Prepare bovine IFN- γ standard by reconstituting contents of vial 4 in 1 ml PBS to give a concentration of 0.5 $\mu\text{g}/\text{ml}$. Leave at room temperature for 15 minutes, then vortex the tube and spin down. Use immediately or store in aliquots at -20°C for future use. We recommend the aliquots not to be refrozen after initial use. For the test, prepare dilutions of the stock using the standard range as a guideline.
 6. Add 100 $\mu\text{l}/\text{well}$ of samples or standards diluted in incubation buffer and incubate for 2 hours at room temperature.
 7. Wash as in step 4.
 8. Add 100 $\mu\text{l}/\text{well}$ of mAb MT307-biotin at 0.25 $\mu\text{g}/\text{ml}$ in incubation buffer. Incubate for 1 hour at room temperature.
 9. Wash as in step 4.
 10. Add 100 $\mu\text{l}/\text{well}$ of Streptavidin-ALP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature.
 11. Wash as in step 4.
 12. Add 100 $\mu\text{l}/\text{well}$ of appropriate substrate solution e.g. p-nitrophenyl-phosphate (pNPP).
 13. Measure the optical density (405 nm for pNPP) in an ELISA reader after suitable developing time.



Distribuito in ITALIA da
Li StarFish S.r.l.
Via Cavour, 35
20063 Cernusco S/N (MI)
telefono 02-92150794
fax 02-92157285
info@listarfish.it
www.listarfish.it